



# DPV: documentation CGM/FGM

Andreas Hungele

Ramona Ranz

As of 20/06/2022 - DPV version 7.50.072

# Contents


<b>1</b>	<b>Documentation in follow-up input mask</b>	<b>2</b>
<b>2</b>	<b>Export of CGM-/FGM-data from the respective software</b>	<b>4</b>
2.1	Abbott – FreeStyle CoPilot . . . . .	4
2.2	Abbott – FreeStyle Libre . . . . .	5
2.3	Abbott - LibreView . . . . .	5
2.4	Dexcom Clarity . . . . .	6
2.5	Dexcom Studio . . . . .	7
2.6	Diasend . . . . .	7
2.7	Glooko . . . . .	11
2.8	Medtronic – CareLink Personal . . . . .	11
2.9	Medtronic – CareLink Pro . . . . .	12
2.10	myDiabby . . . . .	13
2.11	Roche - Accu Chek SmartPix software . . . . .	13
2.12	Tidepool . . . . .	14
<b>3</b>	<b>Import of data in DPV</b>	<b>15</b>
3.1	Import of CGM-/FGM-data . . . . .	15
3.2	Import of basal rate data . . . . .	17
3.3	Import of insulin doses . . . . .	19

# Chapter 1

## Documentation in follow-up input mask

You can find all relevant data for the documentation of CGM-/FGM-data on the tab “metabolic control-self monitoring” in the follow-up input mask.


You can document the amount of days that the patient has used a CGM-/FGM-device in the field “cont. gluc. measurement”. These are the days since the last documentation, maximum for one year.

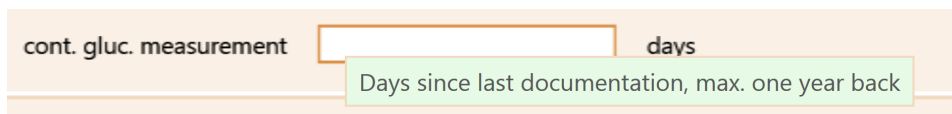

The screenshot shows a form field labeled "cont. gluc. measurement" followed by a text input box and the word "days". A green tooltip box is positioned below the input box, containing the text "Days since last documentation, max. one year back".

Figure 1.1: Input field for amount of days with CGM/FGM


You can copy the amount of days since the last follow-up with the button “cont.” if the patient has used a CGM-/FGM-device without interruptions.

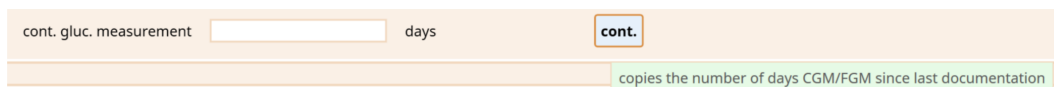

The screenshot shows the "cont. gluc. measurement" field with a "cont." button to its right. A green tooltip box is positioned below the button, containing the text "copies the number of days CGM/FGM since last documentation".

Figure 1.2: Copy amount of days since last follow-up

Please document additional blood glucose measurements in the field “blood glucose / week”.

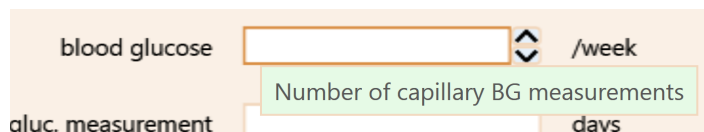

The screenshot shows a form field labeled "blood glucose" followed by a text input box, a dropdown arrow icon, and the text "/week". A green tooltip box is positioned below the input box, containing the text "Number of capillary BG measurements". Below this field, the text "aluc. measurement" and "days" is visible.

Figure 1.3: Number of capillary BG measurements


You can document information about the used CGM-/FGM-device in the region “CGMS”.



The screenshot shows a header with the text 'CGMS' in red. Below it, there are two input fields: 'device' and 'indication', both with dropdown arrows. To the right of the 'device' field is a small icon of a document with a plus sign.

*Figure 1.4: Input fields for device and indication.*

If it is not possible for you to import sensor raw data into DPV, you can enter a few values manually in the region “profiles”. The data relates to the last two weeks.

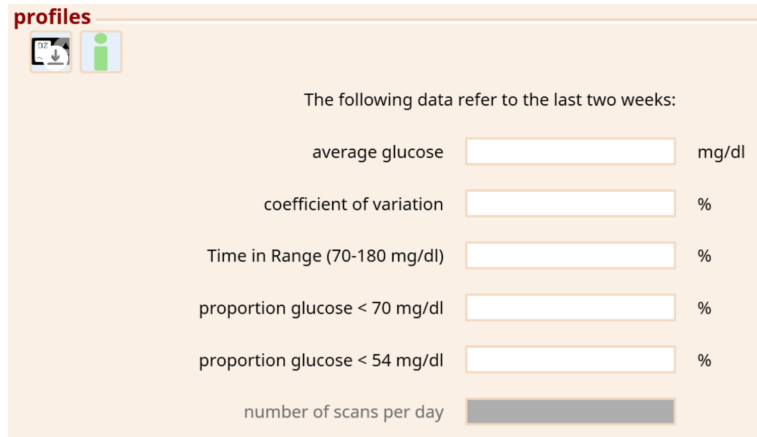

The screenshot shows a section titled 'profiles' in red. Below the title are two icons: a download icon and a person icon. The text 'The following data refer to the last two weeks:' is centered. Below this text are several input fields with labels and units:

Label	Unit
average glucose	mg/dl
coefficient of variation	%
Time in Range (70-180 mg/dl)	%
proportion glucose < 70 mg/dl	%
proportion glucose < 54 mg/dl	%
number of scans per day	

*Figure 1.5: manual documentation of CGM-/FGM-values*

# Chapter 2

## Export of CGM-/FGM-data from the respective software

DPV has no direct access to the CGM-/FGM-devices. Therefore, You have to import the CGM-/FGM-data into the respective software of the CGM-/FGM-manufacturer before you can import the data into DPV.

**Please ensure that the sensor is using the correct date and the correct time before you export the data. Otherwise meaningful analysis of day/night or work day/weekend are not possible.**

**Processing of sensor data saved as a PDF report is not possible.**

Please contact us if one of the following descriptions should not work for you.

### 2.1 Abbott – FreeStyle CoPilot

You can export the measured data via the menu item “File > Export”.

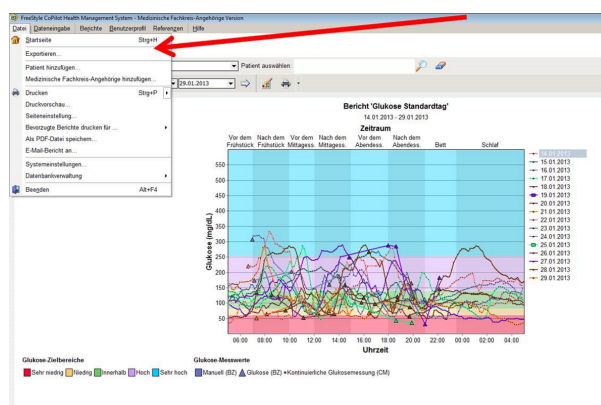


Figure 2.1: Export of the data from FreeStyle CoPilot software

## 2.2 Abbott – FreeStyle Libre

Measured data can be exported from the FreeStyle Libre software into a temporary file via the menu item “File > Export Data”.

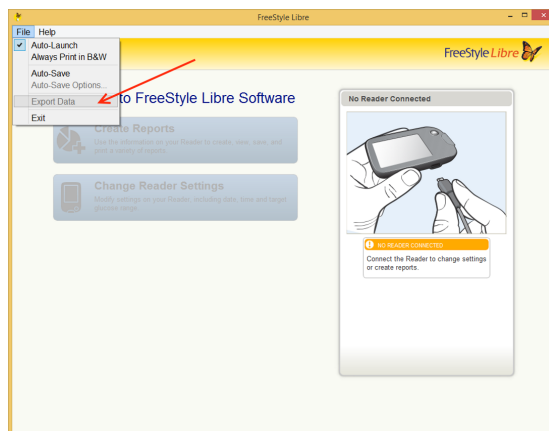


Figure 2.2: Export of the data from FreeStyle Libre software

## 2.3 Abbott - LibreView

After selecting the patient using the patient list, the export of the data can be accessed by clicking the “Profile” link at the top of the left.

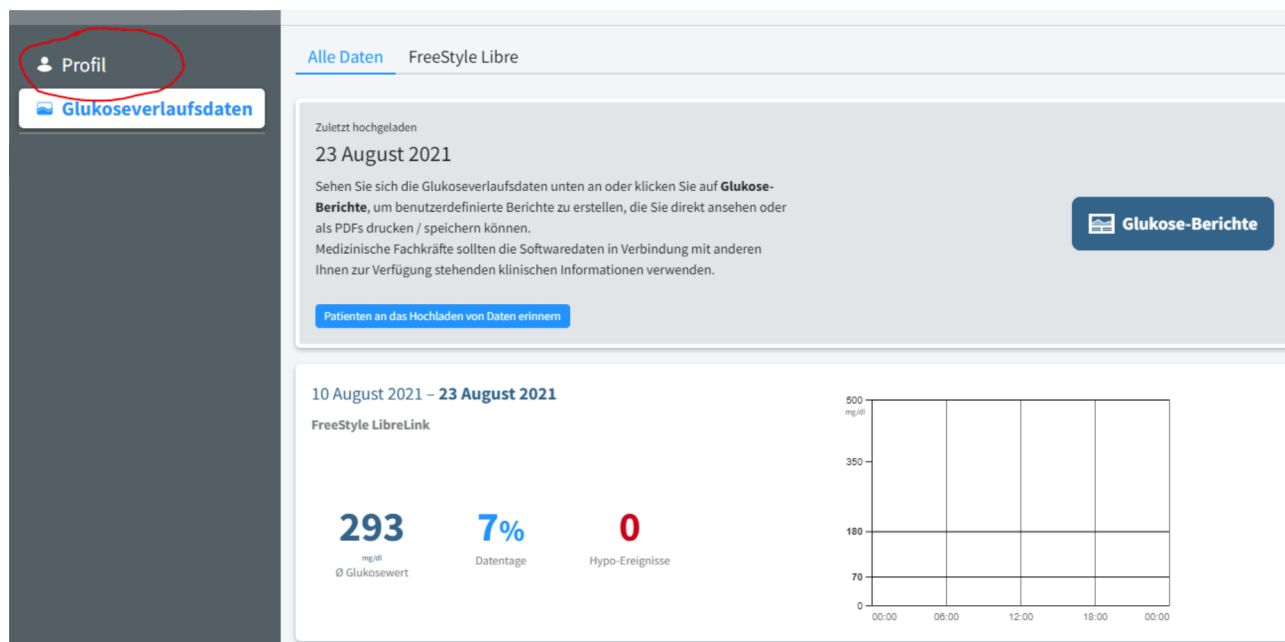


Figure 2.3: Accessing the patient profile

You can download the patient’s glucose values using the link “Download glucose values”.

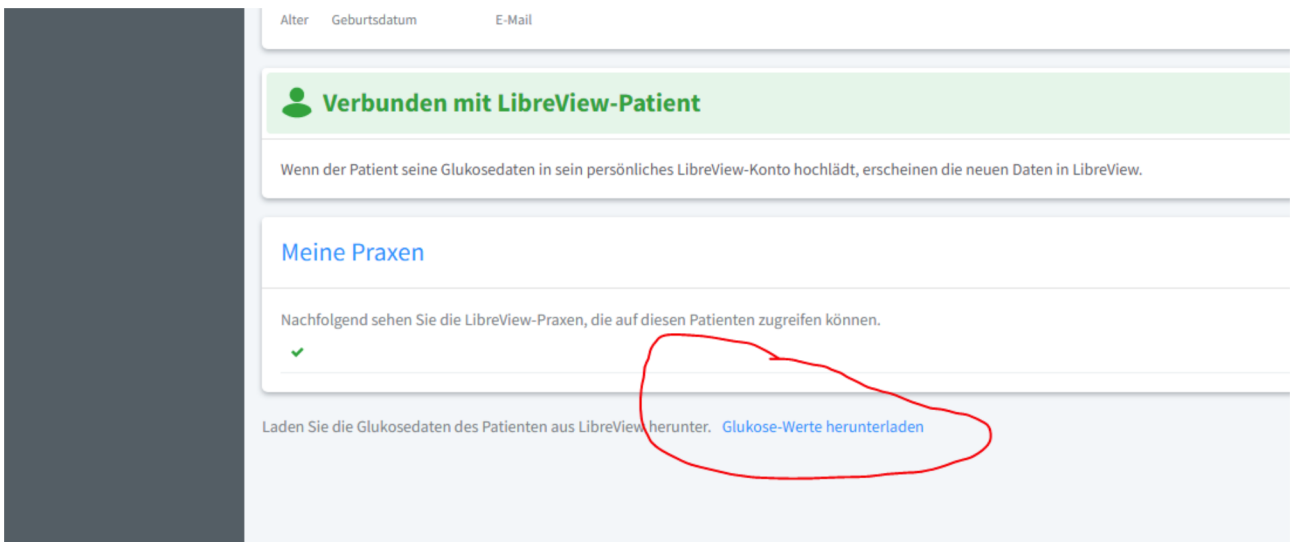


Figure 2.4: Download glucose values

Before the download can be started, you have to verify that you are not a robot. The data can then be downloaded using the "Download" button.

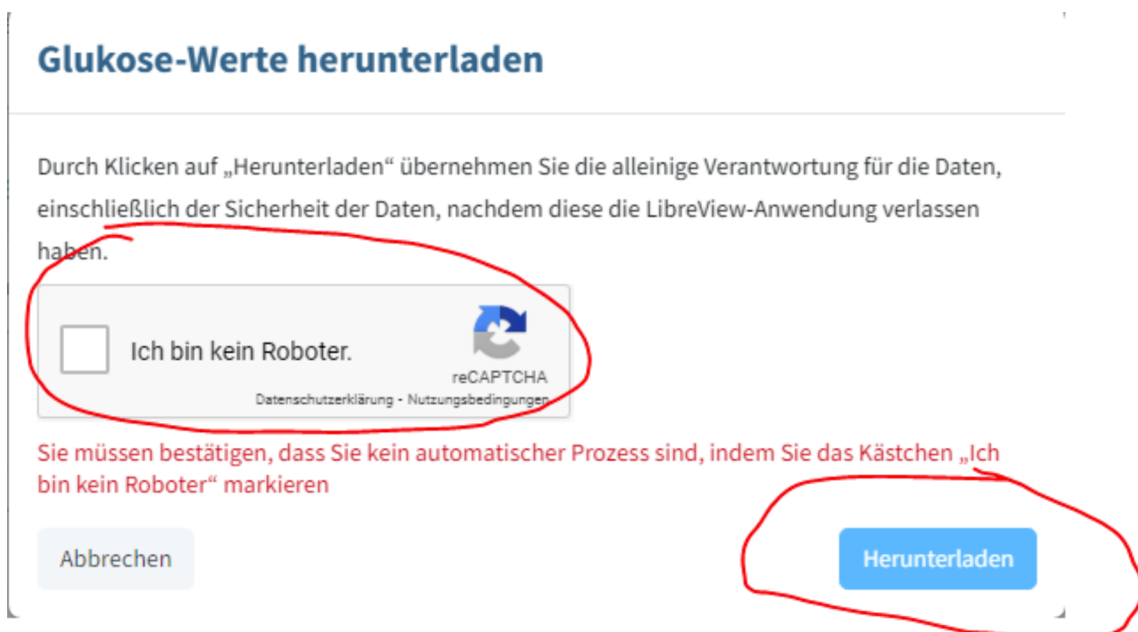


Figure 2.5: Verification and download of glucose values

## 2.4 Dexcom Clarity

You can export the data via "Patient list" > "Export".

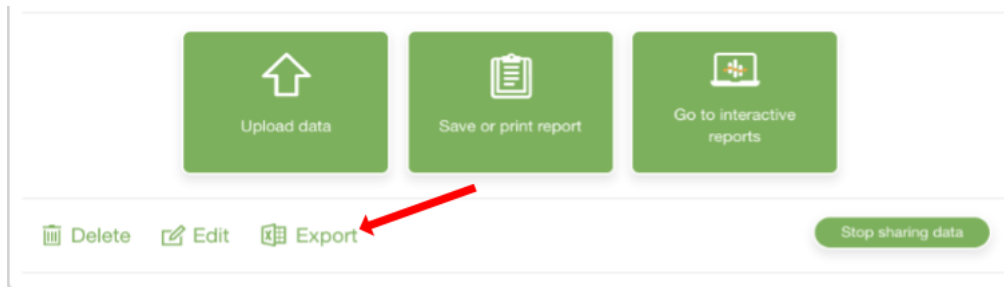


Figure 2.6: Export of the data from Dexcom Clarity

## 2.5 Dexcom Studio

The Dexcom Studio software can be downloaded via <http://www.nintamed.eu/> -> Tab "Service" -> "Downloads" -> "Dexcom G4 Platinum System" -> "Software". You can export the desired patient via the tab item "Patient data" > "Export data". You can choose between the formats "TXT" and "XML". DPV can handle both formats.

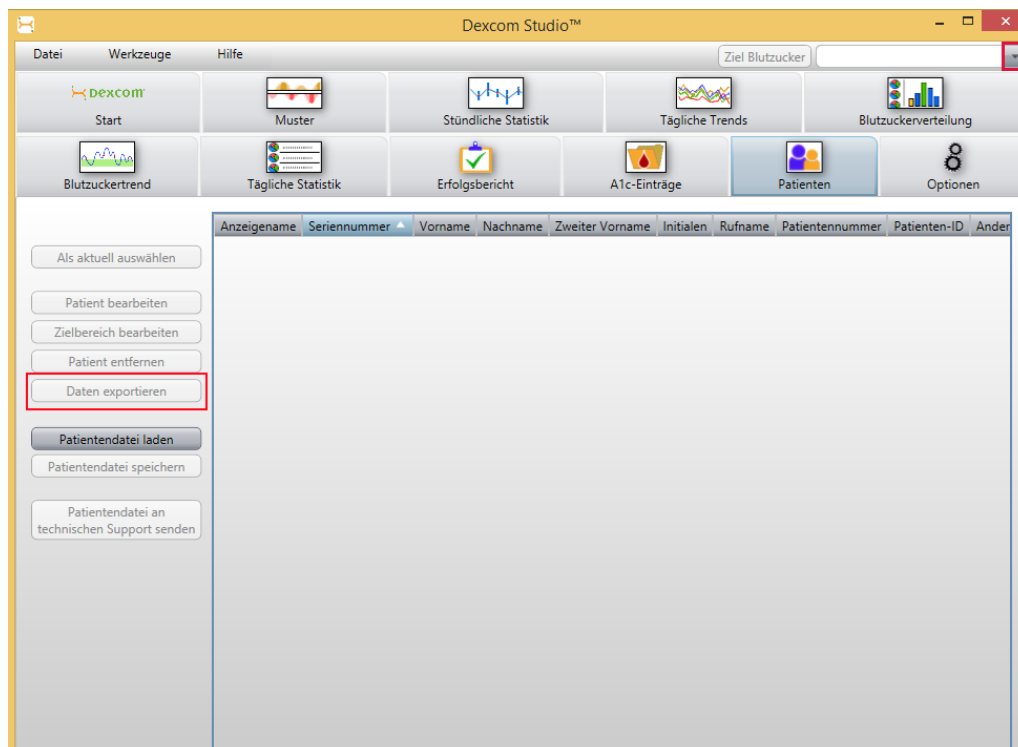


Figure 2.7: Export of the data from Dexcom Studio

## 2.6 Diasend

You can create Excel files for individual patients via "Composition" yourself. Please first select the time interval for which you want to export data. Please then scroll to the bottom of the site. To generate the Excel file please click the button "Export to Microsoft Excel".



November 18, 2015

## Information update: Dexcom G5™ compatibility with Diasend

We understand that some of the previous information regarding the compatibility between Dexcom G5™ and Diasend has caused some confusion. Therefore, here comes a short clarification, the roll-out plan and an FAQ segment. As mentioned earlier, compatibility will be available in pre-determined markets and the method of data transfer varies.

### Data transfer and availability

The new Dexcom G5™ Mobile CGM System has two options. It can be used together with the Dexcom G5™ Mobile Receiver as a “standard CGM unit”. Or it can be used together with an iPhone and the Dexcom app to make the iPhone act as a receiver. This solution is called the Dexcom G5™ Mobile CGM System and the data is wirelessly transferred to the Dexcom cloud.

#### How to upload data from Dexcom G5™ into diasend®

	US	Rest of the world
<b>Cloud to cloud</b> Dexcom G5™ Mobile CGM System - The data is wirelessly transferred via an iPhone and stored in Dexcom’s cloud. This data is then seamlessly integrated into the diasend® cloud by adding the Dexcom G5™ in the connected apps section in your diasend® Personal account.	Available in Q2 2016*	✓
<b>diasend® Transmitter</b> Data from the Dexcom G5™ Mobile Receiver can be uploaded into diasend® by using a micro USB cable and a diasend® Transmitter at a clinic.	✓	✓
<b>diasend® Uploader</b> Data from the Dexcom G5™ Mobile Receiver can be uploaded into diasend® by using a micro USB cable and the diasend® Uploader software.	Within the next month	Within the next month


\* For more information on cloud to cloud availability in the US, please contact Dexcom’s customer service.



# How to upload Dexcom G5™ Mobile - cloud to cloud (Currently not available in the US)


## Requires an iPhone and Dexcom's G5™ Mobile app

If you are using an iPhone to act as a receiver, it is possible to transfer data wirelessly from Dexcom's cloud directly into diasend®.

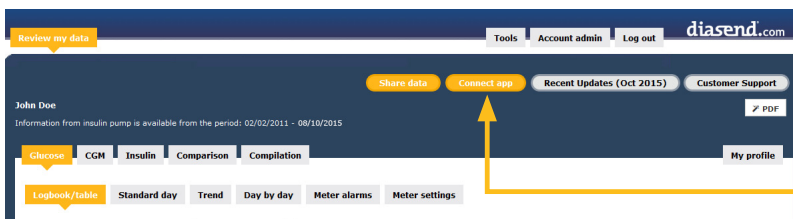

- 1 Log into your diasend® Personal account and click on the button **Connect app**.


**Connect app**

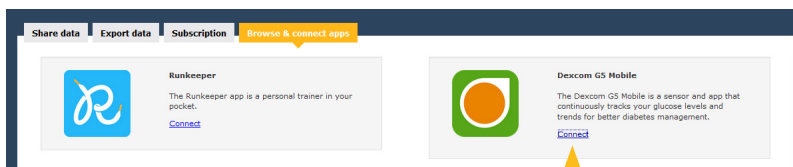

- 2 Click on the **Connect** link for Dexcom G5™ Mobile\*.

**Connect**

\* This feature is currently not available in the US. For more information on the cloud to cloud availability in the US, please contact Dexcom's customer service.

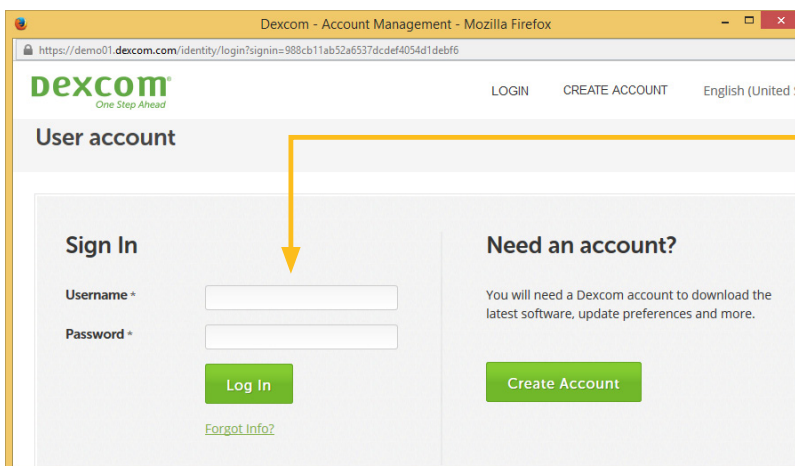

- 3 Sign in with your **Dexcom** username and **password** to authorize and add the Dexcom G5™ Mobile to your diasend® Personal account.



# FAQ

- 1. When can I upload my Dexcom G5™ Mobile Receiver from home (via my diasend® personal account)?**
  - a. The current time frame is to have an update released within the next month.
  - b. In the meantime, if your clinic has a diasend® Clinic account and a diasend® Transmitter, you can upload at the clinic.
- 2. Will my Dexcom G5™ data sync seamlessly with Diasend?**
  - a. Yes, as long as this is available in your country/region and you have an iPhone. The data will sync to diasend® with a 3 hour delay, due to regulatory restrictions.
- 3. I have an Android phone and/or the Dexcom G5™ Mobile System is not yet available in my country. Can I still upload data via diasend®?**
  - a. Yes. Provided that you have a Dexcom G5™ Mobile Receiver and that your clinic has a diasend® Clinic account and a diasend® Transmitter. If so, you will have the ability to upload at the clinic.
- 4. Can my health care provider upload my Dexcom G5™ data?**
  - a. Yes, as long as your health care provider has a diasend® Clinic account and a diasend® Transmitter, then your device can be uploaded at the clinic.
- 5. Can I view my Dexcom G5™ data from home if my receiver was uploaded at the clinic?**
  - a. Yes. As long as you have entered your Dexcom G5™ serial number into your account, under patient profile.
- 6. My Dexcom G5™ data has wirelessly transmitted to my iPhone, but the data does not show up in my diasend® Personal account, why?**
  - a. There is a 3 hour delay for the data to transfer to diasend due to regulatory restrictions.
- 7. I have an iPhone. How do I seamlessly connect my Dexcom G5™ data to Diasend?**
  - a. Login to your diasend® Personal account. Press the **Connect App** button. Then, under the Dexcom G5™ app, press **connect** and enter your Dexcom credentials.

## Any other questions?

Should you have any other questions, please do not hesitate to contact our customer support.

**SUPPORT**  
click here

## 2.7 Glooko

On the “Summary” tab, you can export the patient’s data using the “Export as CSV” link.



Figure 2.8: CSV export via tab “Summary”

In the next step, you can select the period for which the data should be exported. The export of the data can be started using the “Export” button.



Figure 2.9: Selection of period and export

## 2.8 Medtronic – CareLink Personal

You can export the measured data via “Reports > Data export (CSV)”.

- Gesamtüberblick
- Tagesübersicht
- Tagebuch
- Datentabelle
- Standardtag  
Zeitabschnitte
- Standardtag stündlich
- Trendübersicht
- Geräteeinstellungen
- Sensorverläufe Trend-  
Übersicht
- Sensorglukose  
Wöchentlicher Bericht  
nach Zeitabschnitten
- Sensorverläufe  
Übersicht nach  
Mahlzeiten
- Datenexport (CSV)

### Datenexport (CSV)

Startdatum:

Enddatum:

```

1. Last Name: First Name: Date: Time:
2. 20.01.2017 8:00:00 20.01.2017 8:00:00
3. 20.01.2017 8:00:00 20.01.2017 8:00:00
4. 20.01.2017 8:00:00 20.01.2017 8:00:00
5. 20.01.2017 8:00:00 20.01.2017 8:00:00
6. 20.01.2017 8:00:00 20.01.2017 8:00:00
7. 20.01.2017 8:00:00 20.01.2017 8:00:00
8. 20.01.2017 8:00:00 20.01.2017 8:00:00
9. 20.01.2017 8:00:00 20.01.2017 8:00:00
10. 20.01.2017 8:00:00 20.01.2017 8:00:00
11. 20.01.2017 8:00:00 20.01.2017 8:00:00
12. 20.01.2017 8:00:00 20.01.2017 8:00:00
13. 20.01.2017 8:00:00 20.01.2017 8:00:00
14. 20.01.2017 8:00:00 20.01.2017 8:00:00
15. 20.01.2017 8:00:00 20.01.2017 8:00:00
16. 20.01.2017 8:00:00 20.01.2017 8:00:00
17. 20.01.2017 8:00:00 20.01.2017 8:00:00
18. 20.01.2017 8:00:00 20.01.2017 8:00:00
19. 20.01.2017 8:00:00 20.01.2017 8:00:00
20. 20.01.2017 8:00:00 20.01.2017 8:00:00
21. 20.01.2017 8:00:00 20.01.2017 8:00:00
22. 20.01.2017 8:00:00 20.01.2017 8:00:00
23. 20.01.2017 8:00:00 20.01.2017 8:00:00
24. 20.01.2017 8:00:00 20.01.2017 8:00:00
25. 20.01.2017 8:00:00 20.01.2017 8:00:00
26. 20.01.2017 8:00:00 20.01.2017 8:00:00
27. 20.01.2017 8:00:00 20.01.2017 8:00:00
28. 20.01.2017 8:00:00 20.01.2017 8:00:00
29. 20.01.2017 8:00:00 20.01.2017 8:00:00
30. 20.01.2017 8:00:00 20.01.2017 8:00:00
            
```

Der Bericht enthält eine chronologische Liste aller gesammelten Daten im CSV-Format („comma separated value“, durch Komma getrennt).

Dieser Bericht ermöglicht es Ihnen, Ihre Daten zur weiteren Analyse aus dem System zu exportieren.

Berichte erfordern [Adobe® Reader®](#)

Figure 2.10: Export of the data from CareLink Personal

## 2.9 Medtronic – CareLink Pro

You can export the measured data via the button “Export selected data”.

Figure 2.11: Export of the data from CareLink Pro

12

## 2.10 myDiabby

The period for the data export can be selected via “Profile” in the patient data. The period should be set to 3 months. Then the action “CSV export” can be started using “Further actions”.



Figure 2.12: Export of CSV data

In the next window, the type of data to be exported can be selected. Please select the option “Glycaemia, insulins, punctual data” and then click the button “Validate”.

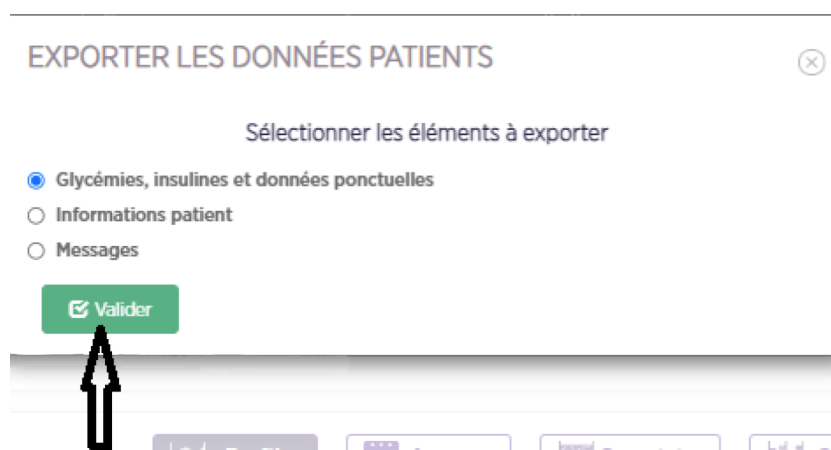


Figure 2.13: Selecting data for export

## 2.11 Roche - Accu Chek SmartPix software

Using the tab “Interfaces” in the system parameters in DPV you can install a plugin for the Accu Chek Smart Pix software which lets you export the sensor data from the Accu Chek SmartPix software. Please note that the Accu Chek SmartPix software must be running during the installation of the plugin.



Figure 2.14: Installing the plugin for Accu Chek SmartPix software

Using the menu item “Device settings > DPVExport” in the patient data you can export the sensor data from the Accu Chek SmartPix software.

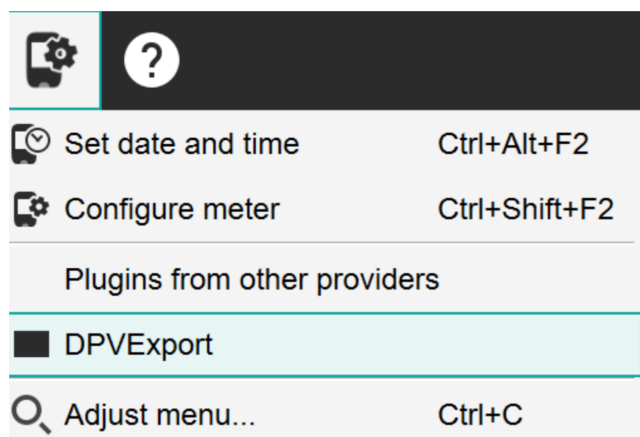


Figure 2.15: Export of the data from the Accu Chek SmartPix software

## 2.12 Tidepool

You can access the settings via the settings button (cogwheel) next to the patient name.



Figure 2.16: Accessing the settings in Tidepool

You can export the data of the patient in the settings in the section “EXPORT MY DATA”. Please select the desired time period and the desired unit for the glucose. The export can be started using the button “Export”.

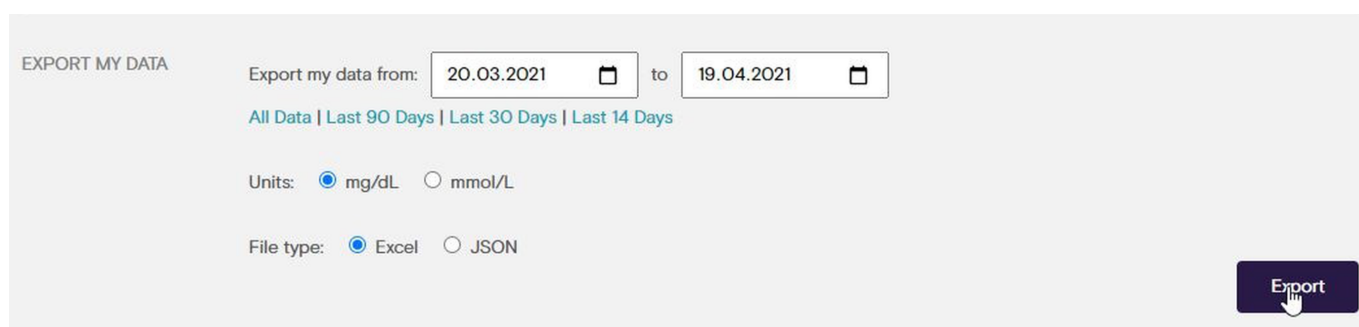


Figure 2.17: Export of the data from Tidepool

Please note that only data exported as “Excel” is currently supported. Data in the format “JSON” can not be processed.

# Chapter 3

## Import of data in DPV

### 3.1 Import of CGM-/FGM-data

You can import adequate files via the button “import CGM-/FGM-data”.



Figure 3.1: Button for the import of CGM-/FGM-data.

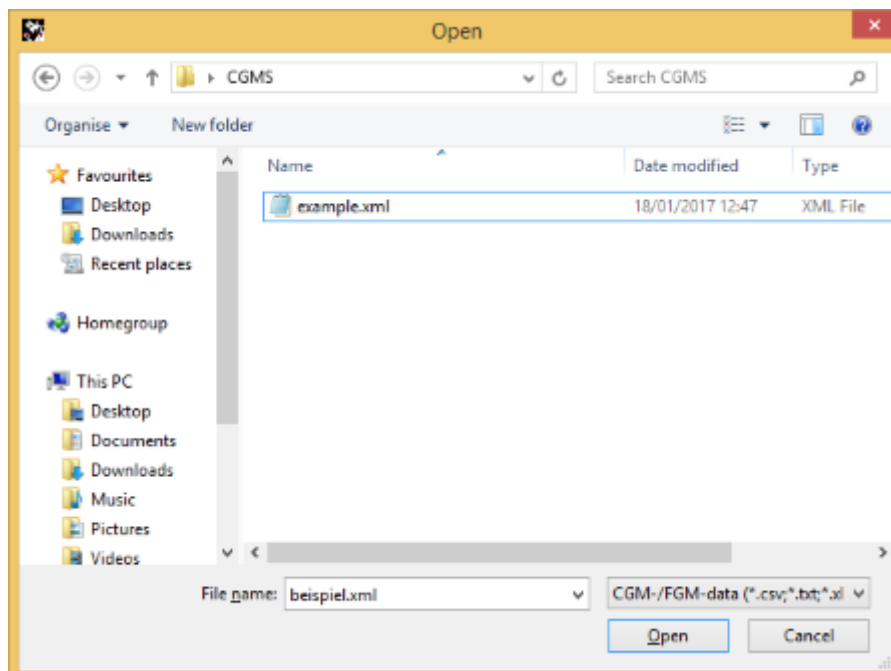


Figure 3.2: Selection of a file with CGM-/FGM-data

After the sensor data has been imported successfully the imported data will be shown in prepared form directly.



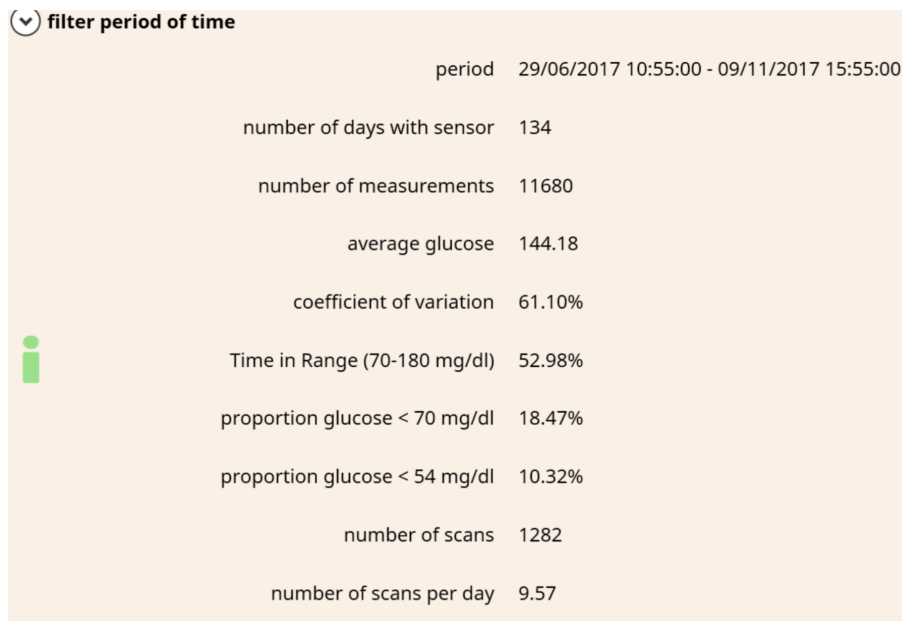


Figure 3.3: Imported data in the CGM input mask

Additionally the sensor data will be displayed in various graphics.

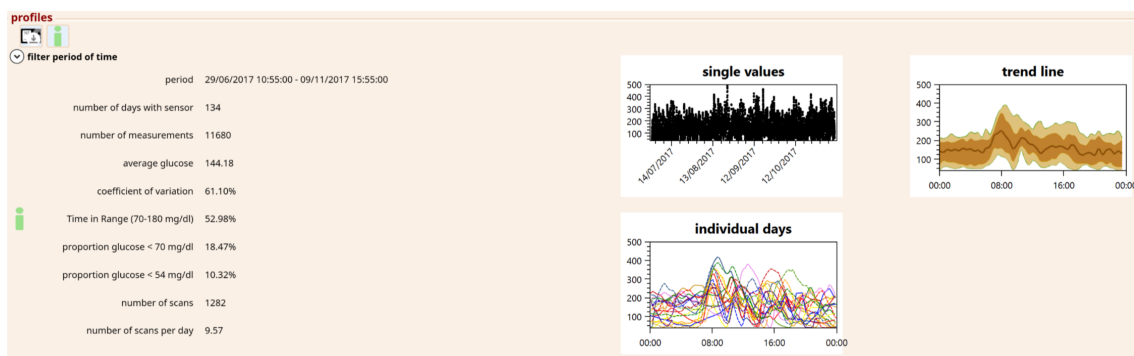


Figure 3.4: Imported data in the CGM input mask with graphics

The following graphics will be displayed:

- **single values:** every measured value will be displayed as a point
- **trend line:** line: median, dark brown surface: 20th - 80th percentile, light brown surface: 5th - 95th percentile, time range: latest 20 days
- **individual days:** layered intraday values of the latest 30 days

Graphics will be displayed enlarged while hovering over the graphic with the mouse.

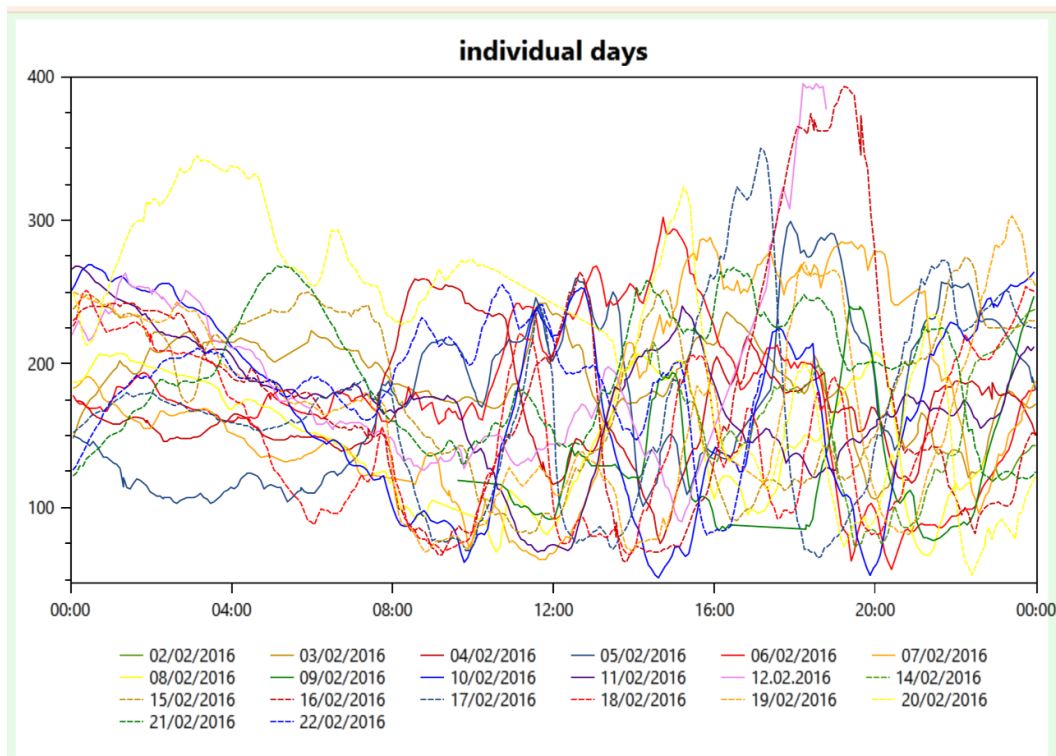


Figure 3.5: Enlarged display of the individual days

Only the most current sensor values related to the current visit date are displayed in the input mask (30 days before examination date). A prompt appears if the imported sensor data does not match the current visit date. The sensor data can then be assigned to another existing visit. Or a new visit can be created directly.


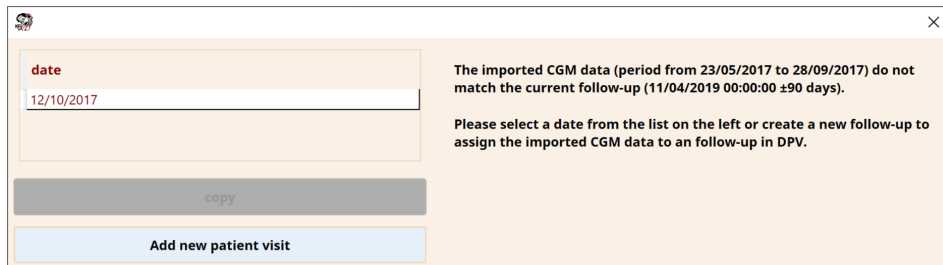


Figure 3.6: Assignment of imported sensor data to another visit

## 3.2 Import of basal rate data

Some of the export files (e.g. CareLink Pro, Diasend, ...) can contain additional data on basal rates. The basal rate data is automatically processed during the import of sensor data.


If basal rate data was processed when importing sensor data, a corresponding message is displayed in DPV.

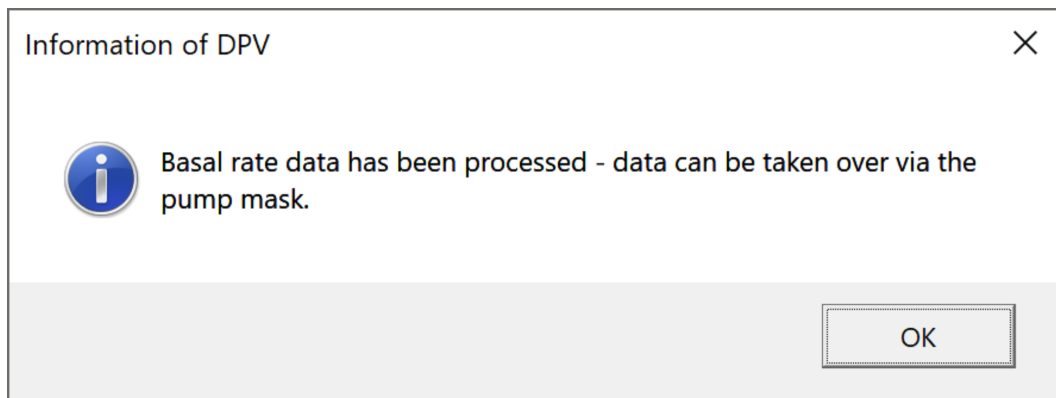

*Figure 3.7: Note, that basal rate data has been processed*

The processed basal rate data can only be imported via the pump input mask. In this way, DPV ensures that the documentation is complete and that there are no inconsistencies in the insulin therapy.

The basal rate data is imported directly when the pump input mask is called. Depending on the type of data imported, the basal rate is imported directly or a selection is displayed for which date the basal rate will be imported.

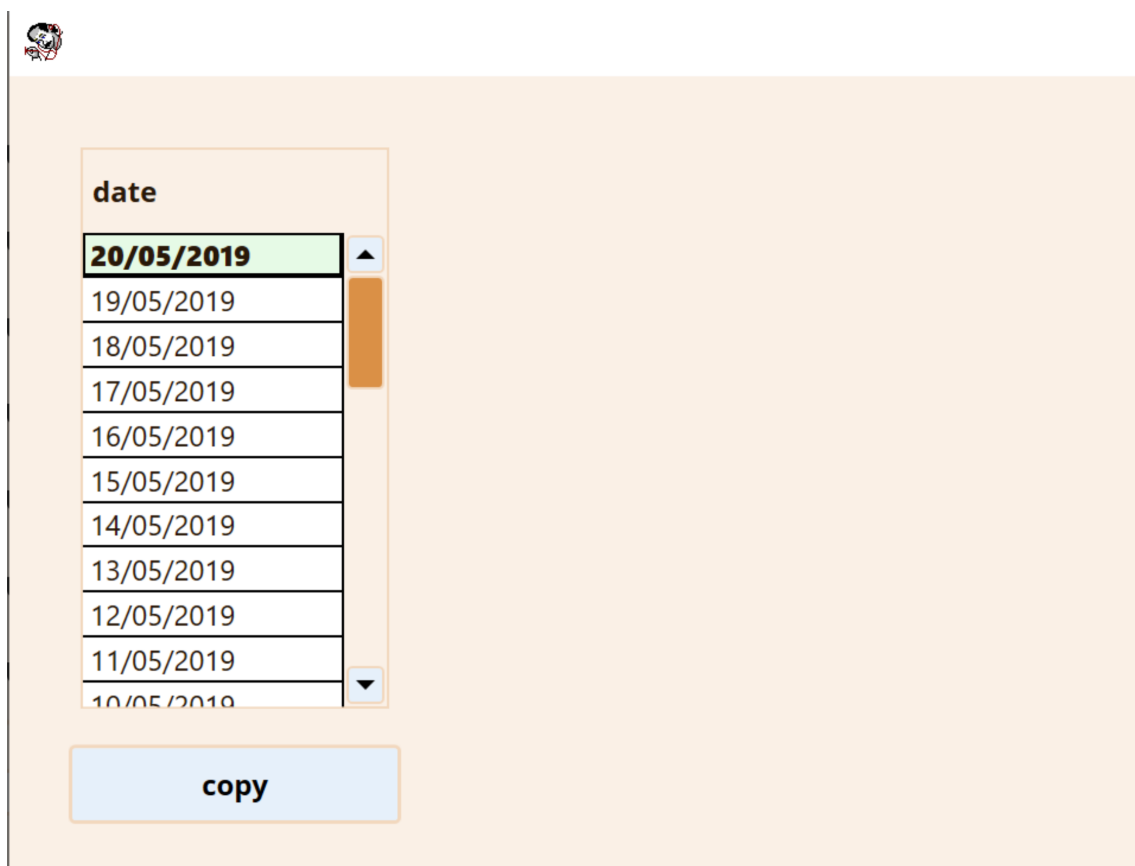


*Figure 3.8: Selection of the date for which basal rate data should be imported*

After selecting the date, you can import the basal rate values by clicking the "copy" button.

### 3.3 Import of insulin doses

Some of the export files (e.g. Glooko) can contain additional data on insulin doses. The insulin doses are automatically processed during the import of sensor data.


If insulin doses have been processed when importing sensor data, a corresponding message is displayed in DPV.

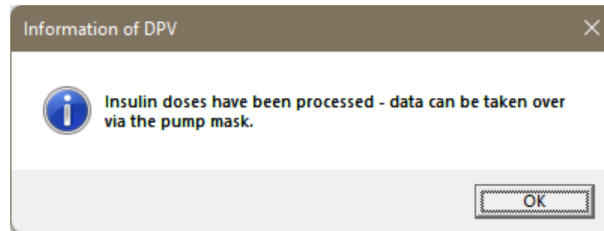

*Figure 3.9: Note, that insulin doses have been processed*

The processed insulin doses can only be imported via the pump input mask. In this way, DPV ensures that the documentation is complete and that there are no inconsistencies in the insulin therapy.

The insulin doses are imported directly when the pump input mask is called. Depending on the type of data imported, the insulin doses are imported directly or a selection is displayed for which date the insulin doses will be imported.

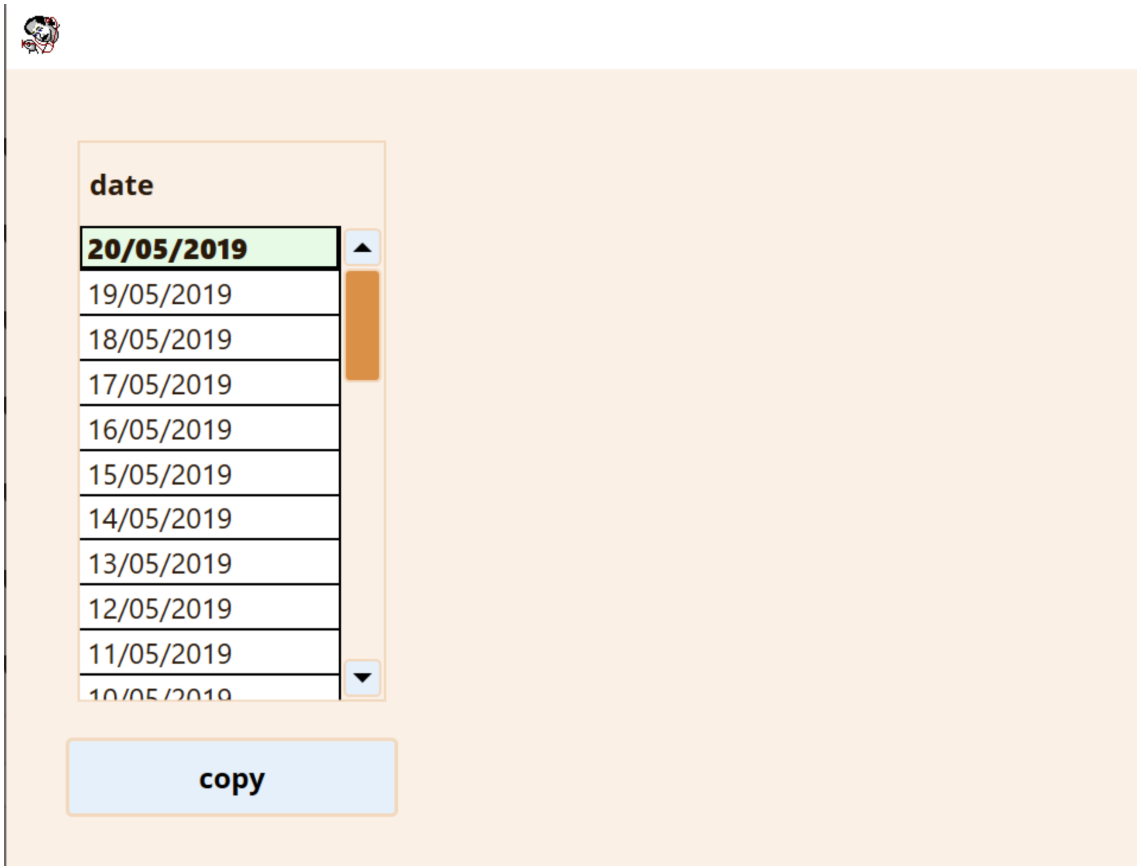

*Figure 3.10: Selection of the date for which insulin doses should be imported*

After selecting the date, you can import the insulin doses by clicking the "copy" button.